




Delivering the Next Generation of “Intelligent” Access Technology.

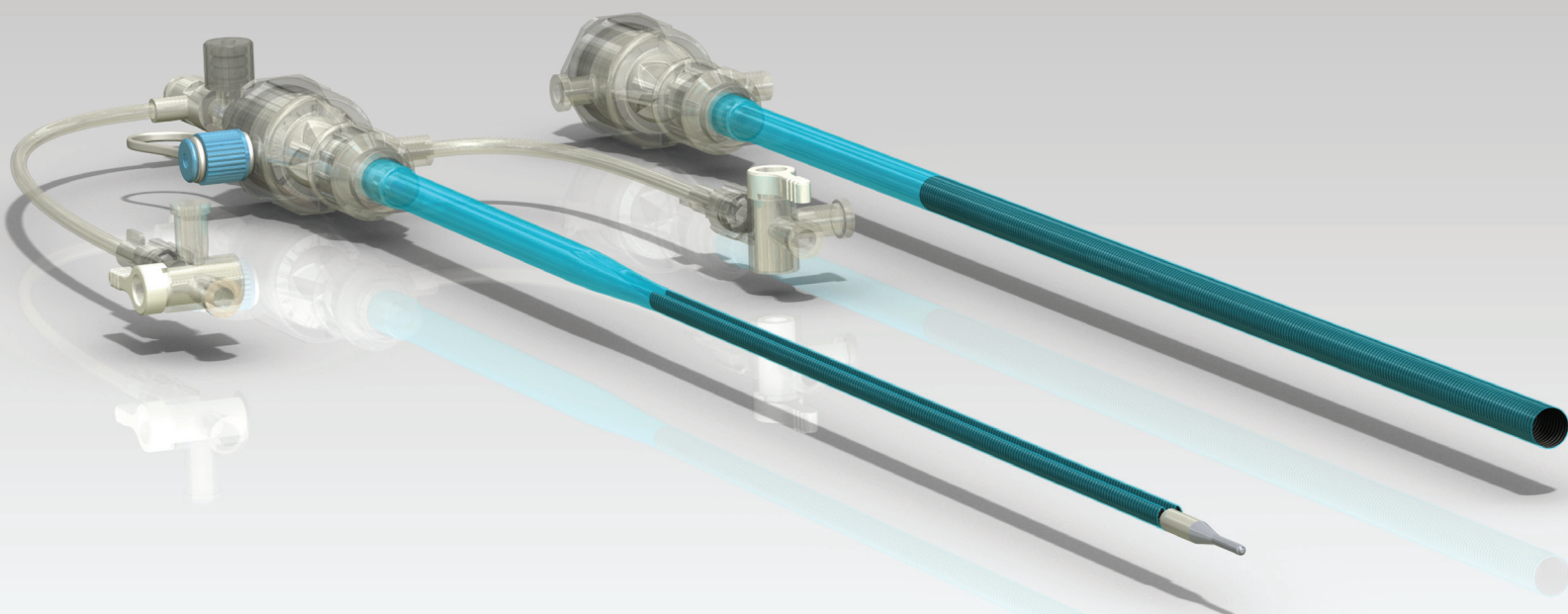

SoloPath™ Balloon Expandable TransFemoral Access System

CDT™ | Controlled Deployment Technology

SoloPath™

SoloPath™ is a Registered Trademark of Onset Medical Corporation

US and International Patents Pending | © Copyright 2010 Onset Medical Corporation | M01003 Rev. C

SoloPath™ Balloon Expandable TransFemoral Access System

CDT™ | Controlled Deployment Technology



Provides the smallest insertion profile of any sheath in its class

- Designed to minimize insertion forces
- Facilitates delivery through difficult anatomy
- Reduced profile enables easier access
- Designed to track vessels similar to a guidewire
- Designed to reduce vessel preparation

Sheath is dilated using Controlled Deployment Technology (CDT™)





- Provides a protective conduit
- Eliminates the use of tapered dilators
- Utilizes natural elasticity of tissue & vessel
- Facilitates more predictable closure

Incorporates state-of-the-art, four-stage hemostasis valve

- Offers maximum hemostasis
- Provides maximum visibility
- Accommodates wire and device
- Offers maximum maneuverability

Unique malleable design

- Designed to protect the vessel
- Provides enhanced maneuverability
- Designed to collapse upon withdrawal

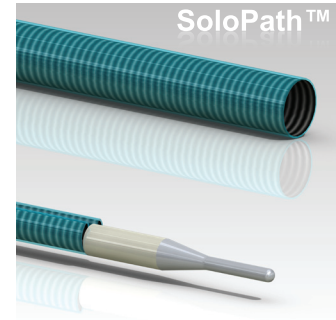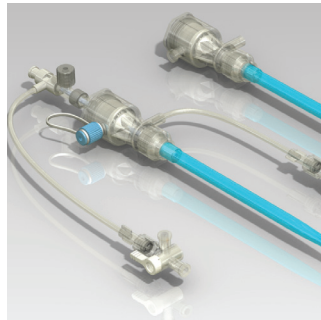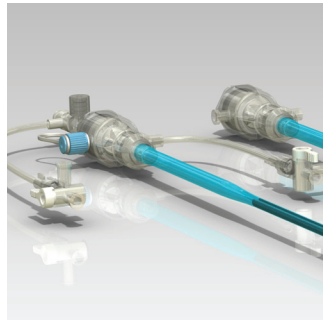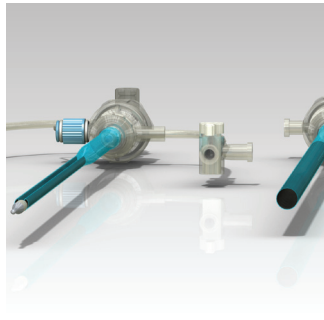

SoloPath™ is a Registered Trademark of Onset Medical Corporation

SOLOPATH™ BALLOON EXPANDABLE TRANSFEMORAL INTRODUCER						
Model	Description	I.D.	O.D.	Expandable Length	Tip Diameter	Working Length
STFI-1825	SoloPath TransFemoral 18Fr x 25cm	18 Fr	21 Fr	20 cm	5.3 Fr	25 cm
STFI-1835	SoloPath TransFemoral 18Fr x 35cm	18 Fr	21 Fr	30 cm	5.3 Fr	35 cm
STFI-1925	SoloPath TransFemoral 19Fr x 25cm	19 Fr	22 Fr	20 cm	5.3 Fr	25 cm
STFI-1935	SoloPath TransFemoral 19Fr x 35cm	19 Fr	22 Fr	30 cm	5.3 Fr	35 cm
STFI-2125	SoloPath TransFemoral 21Fr x 25cm	21 Fr	24 Fr	20 cm	5.3 Fr	25 cm
STFI-2135	SoloPath TransFemoral 21Fr x 35cm	21 Fr	24 Fr	30 cm	5.3 Fr	35 cm
Onset's Part Number	Manufacturer's Part Number	Description			Size	Pressure
C02219	IK1004	Excel Sterile Inflation Device with Stopcock			25cc	30 ATM

Onset Medical Corporation
13900 Alton Parkway Suite 120 | Irvine, CA 92618 USA

Ordering Information
Email: cs@onsetmedical.com | **949.716.1100**